

FAQ's

What are Dermapen Treatments™?

Dp4™ is the breakthrough of the decade utilising the most advanced needle cartridge innovation to vertically pierce the skin. Dp4™ glides over the skin creating millions of fine fractional channels up to 104% faster than other microneedling devices. These channels can carry up to 80% more topical nutrients deeper into the skin. In turn our body's natural response triggers a natural process of repair. The result is an **increased production of collagen and elastin, resurfacing, retexturizing** and overall skin renewal anywhere on the body, from hard to reach areas such as around the eyes, eyelids, earlobes, neck, lips, nose or larger areas of neck and décolleté, stomach and legs for stretch mark treatments.

Dp4™ is the most advanced microneedling device available. Dp4™ treatments are faster, safer, and produce optimal results.

The synergy of Dp4™ and **Dp Dermaceuticals™** used as protocols not only during, but before, and after a Dp4™ microneedling procedures for at-home skincare is what forms a real "Dermapen Treatment.™"

Is the treatment painful?

Traditional dermal rolling techniques proved to be a bloody and painful treatment, requiring the application of topical anaesthetic cream before treatments. Today with Dp4™ there is **no comparison** - this is the **most comfortable**, virtually pain-free microneedling treatment available. Dp4™ procedures are so comfortable there is no need for any therapeutic local anesthesia. Just lay back and enjoy stunning results.

How should I prepare my skin for a Dp4™ procedure?

As the original, trusted leader of microneedling, DermapenWorld™ has collaborated with leading research and development laboratories and clinical experts to create **Dp Dermaceuticals™**, the world's first meso-infusion skincare range exclusively **designed for microneedling procedures**. These active meso-glide combinations are protocols endorsed by thousands of practitioners around the globe for safe infusion into the skin during skin needling procedures.

For optimal results, an appropriate homecare regime of Dp Dermaceuticals™ should be determined by your clinician to encourage the skin to prepare, strengthen and fortify. Continue to use Dp Dermaceuticals™ for a period of 2-3 weeks before and after the Dermapen™ treatments. The Dp Dermaceuticals™ range is suitable for ongoing homecare, a medical grade dermaceutical with powerful "actives" for anti-aging and corrective skincare.

How many treatments are recommended for rejuvenation, acne scarring and stretch marks? (General guide)

Rejuvenation*:	3-6 treatments	6-8 weeks apart
Acne scarring:	4-6 treatments	6-8 weeks apart
Stretch marks:	4-8 treatments	6-8 weeks apart
Pigmentation:	4-6 treatments	2-4 weeks apart

*Superficial treatments maybe recommended at 4 weekly intervals

How soon after filler/muscle relaxant injections can I have a Dp4™ procedure?

While there is no clinical evidence that Dp4™ treatments counteract with any particular injectables, it is important to let your clinician know if you have had an injectable treatment, or are planning one in the near future. **Your clinician** will be able to **advise you** on when to schedule your treatment appropriately.

How long should it take before I see results?

Typically, you will see results after the first treatment. Lasting and more significant results will become visible after 3 to 6 treatments (spaced 2-8 weeks apart). The condition of your **skin** will **continue to improve** over the next 6-12 months after a course of treatments, especially when combined with the recommended post-treatment care.

Is there any downtime after treatments?

Generally no, but you may experience facial redness for 12-48 hours after the treatment. Most people are able to return to normal daily activities immediately after the treatment.

Make an appointment with your
Authorised Treatment Provider today!



Check with your clinician before proceeding with a Dp4 procedure if you are: pregnant, using blood thinners, such as warfarin, or if you are susceptible to keloiding skin. Do not treat over keloid scars, active problematic skin, such as severe solar keratosis, raised moles/warts, eczema or psoriasis. The treatment provider is to be aware this service may have local laws that apply to the operation and use of the Dp4 device. DermapenWorld or it's representatives hold no responsibility for treatment outcomes. Individual results may vary. To be used only as directed.

Dp4™ Dermapen Treatments™ and Dp Dermaceuticals™ are trademarks.
DermapenWorld, Sydney, 2084, Australia. For more information: www.DermapenWorld.com

Dermapenworld™
ORIGINAL • TRUSTED • AUSTRALIAN



**REWRITE
YOUR SKIN'S HISTORY**





With **millions** of treatments performed worldwide you deserve the best

Thousands of practitioners endorse DermapenWorld's world-leading technology which helps you
"Rewrite your skin History"
 without side-effects or downtime.

Doctors say **Dp4**
 can be used to treat the appearance of



Fine Lines & Wrinkles



Acne




Scarring



Sagging Skin




Stretch Marks

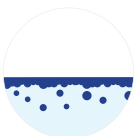

Pigmentation



Enlarged Pores

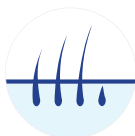

Hair Loss

...and more

Dp4™ & Dp DERMACEUTICALS™

KEY RECOMMENDATIONS FOR OPTIMAL RESULTS



ACNE

Collection to help treat and clear problematic skin and prevent breakouts

In Clinic: Dp4™ with MG-CLR™ & BL™ Mask

Home: CLR Foam Cleanser™, CLR Lotion™ & BL™ Lotion



SCARRING

Effective duo to repair the skin and reduce the appearance of scars

In Clinic: Dp4™ with MG-R.A.S.™

Home: Retinal Active™, SSSCAR™



PIGMENTATION

Collection to address uneven skin tone, discolourations, sun damage.

In Clinic: Dp4™ with MG-BL™ & BL™ Mask

Home: Micro Derm Exfoliant™ & BL™ Lotion



SKIN AGING & SUN DAMAGE

Collection help combat signs of ageing including fine lines, age spots, and skin sagging

In Clinic: Dp4™ with MG-HA35+™ or MG-R.A.S.™ & Hyla Active™ Mask

Home: Antioxidant Cocktail™, Vitamin Rich Repair™ & Retinal Active™



STRETCH MARKS

A powerful combination to reduce the appearance of stretch marks

In Clinic: Dp4™ with Hyla Active™, Retinal Active™ and SSSCAR™

Home: Hyla Active™, Retinal Active™ & SSSCAR™



Ask your treatment provider about the full range of Dp Dermaceuticals™.

Your Treatment Provider will assess your skin and treatment area to recommend a treatment plan suitable for you.